

Rabbit anti-Human Cholecystokinin Type A Receptor (CCK1) (AA 409-428) (CT)

#Cat: NB-49-003-50ul

Size: 50ul


Product Description: CCKA antibody, human, C-term.
CCK1
Cholecystokinin type A receptor
Cholecystokinin type 1 receptor
G Protein-coupled receptor
Polyclonal antibody raised against the
C-terminus of the human CCKA receptor

Immunogen: residues 409-428 coupled to KLH (NP_000721)
KLH- tgaslrfsysyhmsasvppq

Host: Rabbit

Form: serum with 0.02% NaN₃ as preservative

Quantity: 50 µl

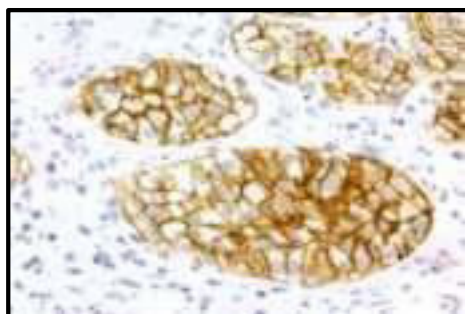

CCKA receptor immunoreactivity in a human carcinoid tumor.

Applications: ICH (p), IC, W
W = Western blot (1:5000)
IC = Immunocytochemistry (1:5000)
ICH (p) = Immunohistochemistry of formalin-fixed, paraffin-embedded tissues (1:2500)

Species cross-reactivity: human

Sensitivity: 50 ng

Storage 4 °C short term
-20 °C long term